

Inspiring for Splendid Healthcare

ISO
9001:2008

ISO
13485:2012



PVA DRESSINGS

For EAR & NOSE

OSSEOUS

PVA DRESSINGS



EON

Inspiring for Splendid Healthcare

MEDITECH

PVT. LTD.




Ear PVA Dressings

Ear Wick:

Osseous ear wicks are designed to aid the treatment of otitis externa. They deliver medication, such as an antibiotics solution, to the ear canal and the ear canal membrane

- ✦ Osseous ear wicks are designed to deliver antibiotics during the treatment of otitis externa
- ✦ Perfect Design to use in treatments for otitis externa
- ✦ Delivers medication to ear canal
- ✦ Fenestrated variety available
- ✦ Soft and resilient

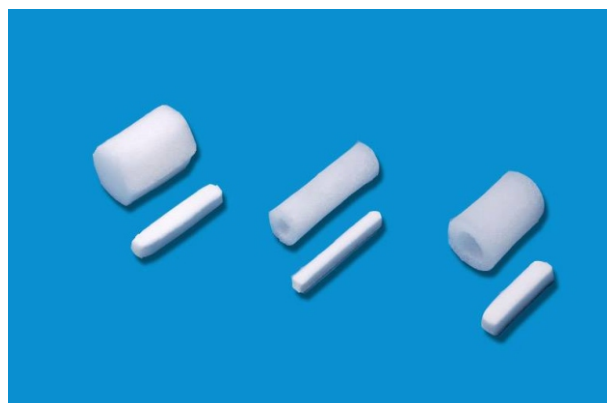


REF Number	Product Name	Specification	Product Description	No. Of Piece/Package
E7161	Ear Wick	Dia.9.0mmX15.0mm	No String	20 Pcs, Gamma Sterile
E7162	Ear Wick	Dia.9.0mmX24.0mm	No String	20 Pcs, Gamma Sterile
E7163	Ear Wick	Dia.12.0mmX15.0mm	No String	20 Pcs, Gamma Sterile
E7164	Ear Wick	Dia.12.0mmX24.0mm	No String	20 Pcs, Gamma Sterile
E7169	Ear Wick	Dia.9.0mmX15.0mm	Fenestrated	20 Pcs, Gamma Sterile
E7170	Ear Wick	Dia.9.0mmX24.0mm	Fenestrated	20 Pcs, Gamma Sterile



Ear PVA Dressings

Ear Pack:

- ✦ Designed for use in otic surgery, Osseous ear packs help absorb fluids during surgery. The soft and porous PVA foam material expands evenly to support the ear canal after surgery. Osseous ear packs are sterile and come with rounded edges for quick and easy insertion
- ✦ Osseous ear packs are designed for use in otic surgery. They absorb fluids and support canals after surgery
- ✦ Helps prevent stenosis of the ear canal
- ✦ Porous foam for rapid absorption

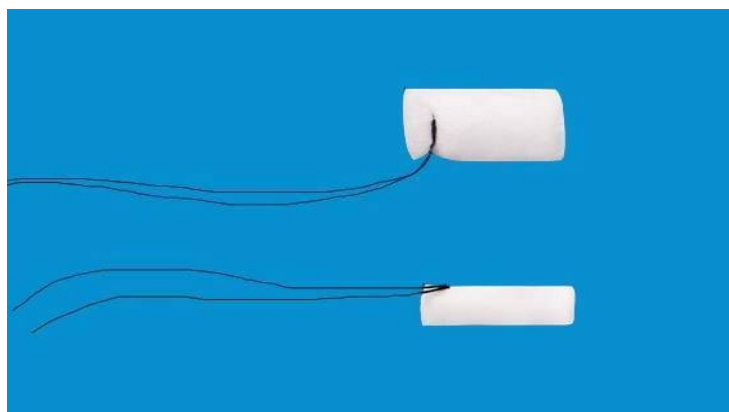


REF Number	Product Name	Specification	Product Description	No. Of Piece/Package
E7165	Ear Pack	Dia.9.0mmX15.0mm	With String	20 Pcs, Gamma Sterile
E7166	Ear Pack	Dia.9.0mmX24.0mm	With String	20 Pcs, Gamma Sterile
E7167	Ear Pack	Dia.12.0mmX15.0mm	With String	20 Pcs, Gamma Sterile
E7168	Ear Pack	Dia.12.0mmX24.0mm	With String	20 Pcs, Gamma Sterile



Standard Nasal Dressings

Osseous standard nasal packs are designed for use in a variety of nasal surgeries: septoplasty, turbinectomy, rhinoplasty, etc. They fill and support nasal cavities during surgery, allowing smooth septum ventilation to ensure patient comfort. Osseous nasal packs also absorb post-operative bleeding effectively and efficiently

- ✦ Designed for nasal surgeries.
- ✦ Gentle compression to stop bleeding and absorb blood
- ✦ Excellent alternative to cumbersome gauze packing

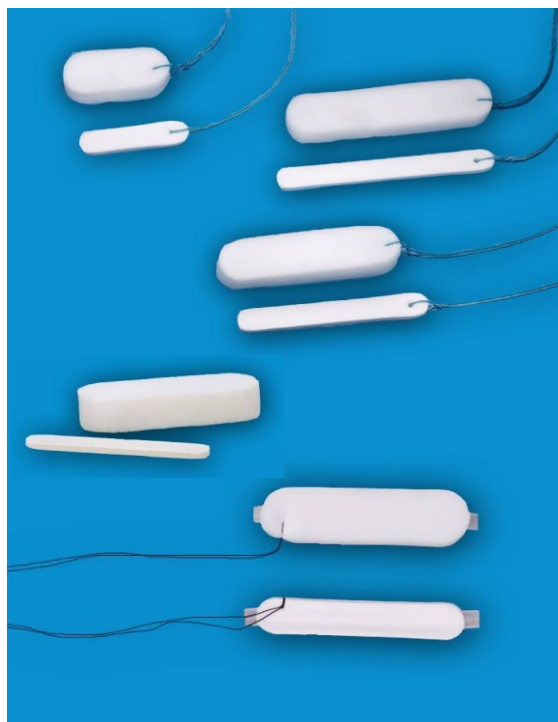


REF Number	Product Name	Specification	Product Description	No. Of Piece/Package
E7176	Nasal Dressing	4.5cmX1.5cmX2.0cm	No String	10 Pcs, Gamma Sterile
E7177	Nasal Dressing	4.5cmX1.5cmX2.0cm	With String	10 Pcs, Gamma Sterile
E7178	Nasal Dressing	4.5cmX1.5cmX2.0cm	With String & Tube	10 Pcs, Gamma Sterile
E7179	Nasal Dressing	4.5cmX1.0cmX2.0cm	Thin	10 Pcs, Gamma Sterile
E7180	Nasal Dressing	8.0cmX1.5cmX2.0cm	No String	10 Pcs, Gamma Sterile
E7181	Nasal Dressing	8.0cmX1.5cmX2.0cm	With String	10 Pcs, Gamma Sterile
E7182	Nasal Dressing	8.0cmX1.5cmX2.0cm	With String & Tube	10 Pcs, Gamma Sterile
E7183	Nasal Dressing	8.0cmX1.0cmX2.0cm	Thin	10 Pcs, Gamma Sterile



Epistaxis Nasal Dressings

Osseous epistaxis nasal packs are used in the treatment of anterior and posterior nose bleeding. The soft round edges create less friction and reduce irritation during nasal passage insertion. The epistaxis packs expand to gently compress and effectively stop nose bleeding

- ✦ Epistaxis packs stop nosebleeds, safely and effectively
- ✦ Designed to treat nosebleeds
- ✦ Softness reduces irritation
- ✦ Quick and uniform expansion
- ✦ Fluid absorption causes the packing to expand producing a tamponading effect

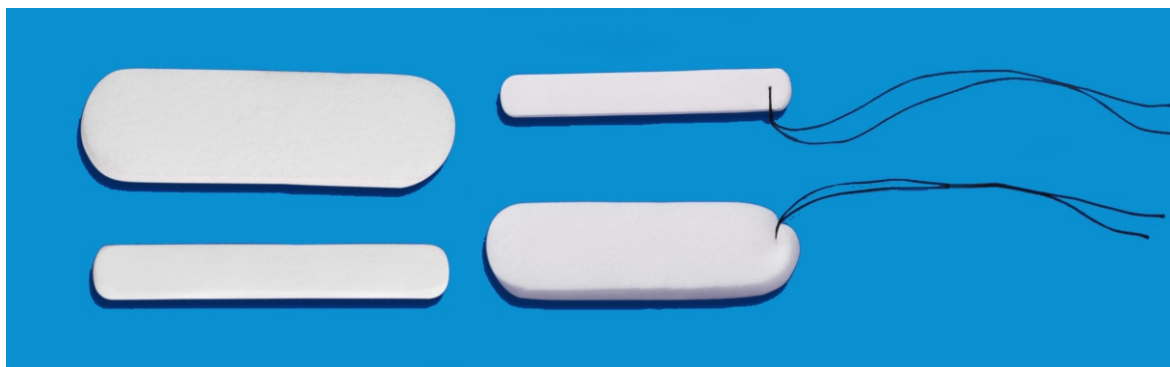


REF Number	Product Name	Specification	Product Description	No. Of Piece/Package
E7184	Epistaxis Nasal Dressing	10.0cmX1.5cmX2.5cm	No String	10 Pcs, Gamma Sterile
E7185	Epistaxis Nasal Dressing	10.0cmX1.5cmX2.5cm	With String	10 Pcs, Gamma Sterile
E7186	Epistaxis Nasal Dressing	10.0cmX1.5cmX2.5cm	With String & Tube	10 Pcs, Gamma Sterile
E7187	Epistaxis Nasal Dressing	10.0cmX1.0cmX2.5cm	Thin	10 Pcs, Gamma Sterile
E7193	Epistaxis Nasal Dressing	5.5cmX1.0cmX2.5cm	Pediatric	10 Pcs, Gamma Sterile



Anatomical Nasal Dressings

Designed to fill the nasal cavity, Osseous anatomical nasal packs are used during and after nasal surgery. With optional airway tube, they enhance patient comfort by ensuring easy breathing. Moderately-sized, and with soft yet supportive texture for easy placement and removal, Osseous anatomical packs effectively reduce patient discomfort

- ✦ Anatomical packs are used during and after nasal surgery to support nasal cavities and facilitate proper breathing
- ✦ Support the nasal cavity after surgery
- ✦ Excellent alternative to cumbersome gauze packing

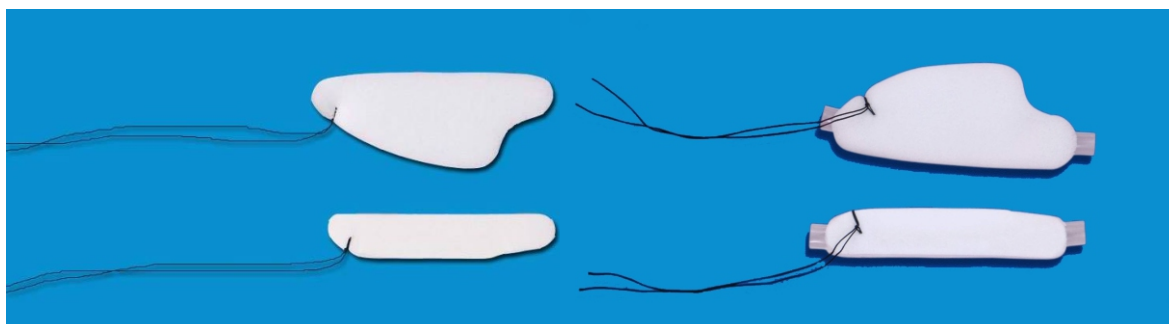


REF Number	Product Name	Specification	Product Description	No. Of Piece/Package
E7188	Anatomical Nasal Dressing	8.0cmX1.5cmX3.0cm	With String	10 Pcs, Gamma Sterile
E7189	Anatomical Nasal Dressing	8.0cmX1.5cmX3.0cm	With String & Tube	10 Pcs, Gamma Sterile
E7190	Anatomical Nasal Dressing	8.0cmX1.5cmX2.0cm	Thin	10 Pcs, Gamma Sterile



Sinus Packs

Osseous sinus packs are made with soft, porous, resilient PVA foam. They are designed to support the sinus cavity and the osteomeatal complex to ensure airflow through nose. They also expand to stop bleeding through gentle compression, and facilitate faster recovery. They have rounded edges to allow easy insertion during surgery, and provide comfort even during removal.

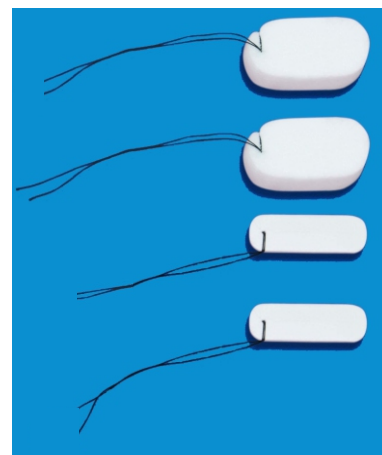

- ✦ Sinus packs are designed for use after sinus surgery, to ensure proper breathing.
- ✦ Ideal for FESS (Functional Endoscopic Sinus Surgery).
- ✦ Reduce bleeding and facilitate healing.
- ✦ Allows sinus cavity to remain open while absorbing post-operative fluids.
- ✦ Fine pore construction reduces tissue adherence.

REF Number	Product Name	Specification	Product Description	No. Of Piece/Package
E7191	Sinus Pack	3.5cmX1.2cmX0.6cm	With String	10 Pcs, Gamma Sterile
E7192	Sinus Pack	3.5cmX1.2cmX0.9cm	With String	10 Pcs, Gamma Sterile

Common Advantages for all kind of Osseous PVA dressings:

- ✦ Lint and fiber-free
- ✦ Extremely hydrophilic
- ✦ Soft and comfortable for patient use
- ✦ Excellent tensile strength, Controllable range of pore size
- ✦ Packaged in clean-room facilities
- ✦ Excellent alternative to cumbersome gauze packing
- ✦ All products have good biocompatibility, low antigenicity, tissue reaction and Transmissibility